

## MAG 1000 HIGH FREQUENCY MAGNETOTHERAPY - 20 PROGRAMMES

CE 0068



28338

### • 28338 MAG 1000 MAGNETOTHERAPY - 2 channels - 20 programmes

Mag 1000 is a high frequency PEMF device particularly recommended for the treatment of pathologies affecting soft tissues like muscles and tendons. It offers 20 pre-set programmes with 3 adjustable intensity levels.

The device generates an electric signal with a carrier frequency of 27 MHz, which is modulated with specific frequencies (from 8 Hz to 640 Hz) and therefore is able to soothe the anatomical site affected by pain. It is projected and indicated for healing, rehabilitation and functional recovery treatments of many pathologies. Continuous operating mode. It is intended for both professional and homecare use.

Size: 153x91x38 mm

Supplied with elastic therapeutic belt, carrying bag, rechargeable battery pack 2,300mAh, battery charger, user manual in GB, FR, IT, ES, DE. Medical Device Class IIa.



28338

### 20 PROGRAMMES

Arthritis-Arthrosis	Bruises
Rheumatisms	Myalgia
Pain-articular prosthesis	Tendonitis
Herniated disc	Hematoma
Muscular pain	Epicondylitis
Osteoporosis	Epytrocleitis
Fractures	Lumbago
Sprains	Muscular contractures
Dislocations	2h Medium application
Traumas	4h Long application

## TAP 2000 MAGNETOTHERAPY MAT

### • 28339 TAP 2000 MAGNETOTHERAPY MAT for MAG 2000

TAP 2000 is a 40x40 cm mat with 2 cables of 3 solenoids each (total 6 solenoids), ideal for magnetotherapy treatments.


It is compatible with high intensity/low frequency pulsed electromagnetic field devices: MAG 2000, MAG 2000 Plus and MAG 2000 Premium.

Supplied with user manual in GB, FR, IT, ES, DE.

Class I Medical Device.

CE

28339

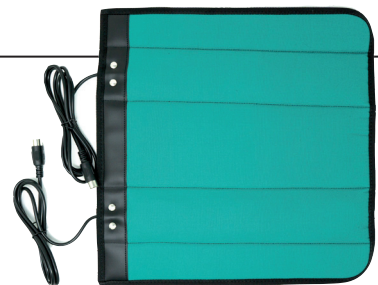

## KURA TENS - 2 CHANNELS

### • 28431 KURA TENS - 6 programmes - GB, FR, ES, PT, DE, NL, PL, CZ, SE, Arabic

### • 28432 KURA TENS - 6 programmes - IT, DK, LT, FI, HU, RO, HR, SI, GR, BG

Dual channel electrotherapeutic device with massage feature ideal for pain relief associated with aching muscles due to physical exertion or regular work and domestic activities.

It is also aimed at providing symptomatic relief and managing chronic and intractable pain, as well as pain relief related to arthritis. Medical Device Class II a.

- 2 channels with 4 electrode pads

- 15 intensities per channel


- 6 main modes (composite, thump, press, slap, acupuncture, muscular relaxation)

- 24 sub-modes (4 for each main mode)

- large LCD with blue backlight.

Supplied with 2 wires (120 cm), 4 electrode pads (50x50 mm) and multilingual manual (see above).

### • 28408 REUSABLE GELLED ELECTRODES 50x50 mm - bag of 4

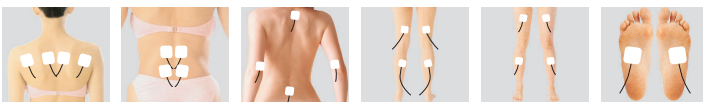

28431, 28432

CE 0197



### TECHNICAL SPECIFICATIONS

Output channels:	2
Body parts:	waist, back, arm, leg, joint, shoulder, adductor and gluteus, chest, belly, feet
Treatment duration:	15 to 30 minutes
Band width:	0~100 us
Pulse frequency range:	0-100 Hz
Output voltage:	0-60 V ± 20% (500 ohm load)
Power supply:	4 x 1.5 V AAA batteries
Size - Weight:	120.3x60.3x20.6 mm - 73 g (batteries excluded)
Classification:	internally powered, Type BF applied part, IP22

## MADRID SPHYGMOMANOMETER

### • 50000 MADRID ANEROID SPHYGMOMANOMETER

Professional, high quality sphygmomanometer with 0 to 300 mmHg scale. Two handling modes for right and left-hand users.

Large dial Ø 59.5 mm with PVC shockproof ring.

Zinc-alloy handle with nickel plated. Copper male and female connector.

Supplied with adult nylon cuff (22-32 cm) with D-ring, 1-tube PVC bladder, ribbed PVC bulb, black nylon bag and multilingual box and manual.



Luminescent dial



Male and Female connector



Left and Right hand use



50000

CE 0123